

# DuraClip

## Installation Manual




**U.S.:** 700 Ensminger Rd., Suite 112, Tonawanda, NY 14150  
TOLL FREE TEL: 1.888.433.6818 TOLL FREE FAX: 1.888.433.6834

**Canada:** 84 Citation Dr., Concord, ON L4K 3C1 Canada  
TOLL FREE TEL: 1.866.748.7943 FAX: 905.695.0138



# DURA CLIP

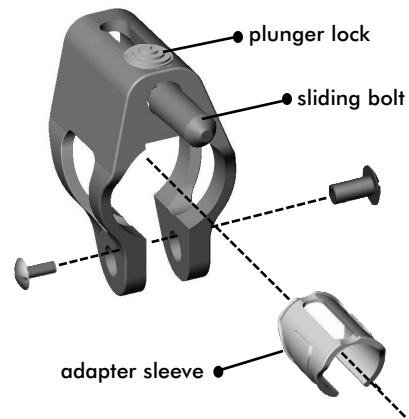

The **DuraClip** is a “lock and release” device designed to secure the drop hooks from a wheelchair (drop-base) cushion onto the wheelchair seat frame/tube. By engaging the plunger lock mechanism, the sliding bolt holds the drop hook (seat) in place; By retracting the plunger lock, the seat can be easily removed for wheelchair transportation.

## INSTALLATION INSTRUCTIONS:

### Standard Application (using two DuraClips):

1. Position the wheelchair seat in the desired position on the seat frame tubes and snap the seat (drop hooks) into place.

2. Working at the front of the wheelchair, install a DuraClip onto each seat tube (left & right) as illustrated in image A.

*\*note: the DuraClips must be installed in front of the drop hooks such that the sliding bolt extends rearward over top of the drop hooks- see image C.*

3. If the wheelchair uses 7/8" (22mm) diameter seat tubes, install adapter sleeves as indicated in image B.

4. Slide each DuraClip along the seat tube until it rests up against the drop hook. Install the mounting hardware provided, and secure<sup>2</sup> the DuraClip in position on the seat frame/tube-see image C.

*\*note: tighten hardware until clip is snug on the seat tube- **do not over-tighten!***

5. With both DuraClips installed, press down on the plunger lock mechanism and extend the sliding bolt over top of the drop hooks to secure the seat on the wheelchair frame. (Ensure the plunger mechanism is locked into position).

6. To unlock/release the seat from the wheelchair, depress the plunger mechanism and retract the sliding bolt on each DuraClip. Remove seat (drop base) cushion.

