



HEAD SUPPORT OPTIONS:

STANDARD HARDWARE

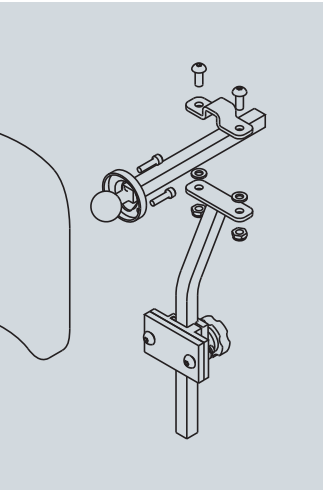


- > Comes standard with Otto Bock, Standard & Stealth Comfort Plus Pads

STEALTH "WORLD'S BEST"

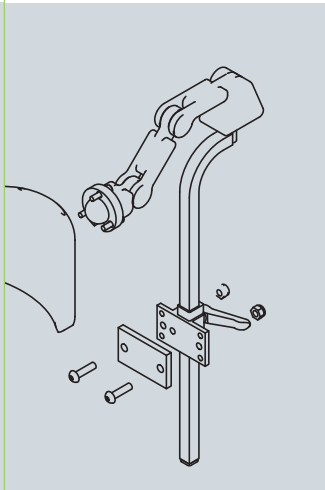


- > Upgrade option for Otto Bock, Standard & Stealth Comfort Plus Pads
- > Dual pivot balls for pad orientation
- > Depth adjustable through dual link system
- > Height adjustable
- > Includes mounting clamp for attachment to back with latch style release for headrest removal

MULTI-AXIS

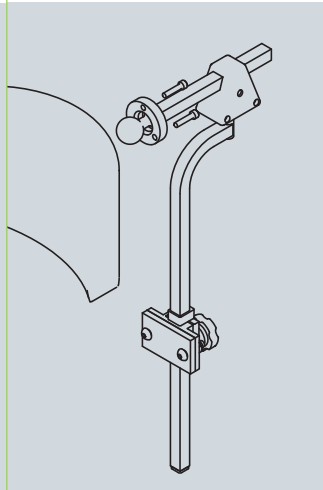

- > Upgrade option for Otto Bock, Standard & Stealth Comfort Plus pads
- > Dual pivot balls for pad orientation
- > Depth adjustable through horizontal slide bar
- > Height adjustable
- > Includes mounting clamp for attachment to back with knob style release for headrest removal

FIXED MOUNTING HARDWARE

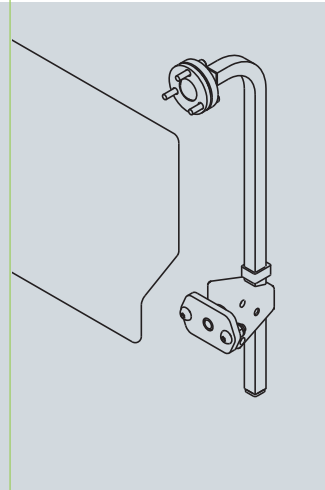

- > Dual pivot ball for pad orientation
- > Limited adjustability – recommended for Contour or High backs
- > Height adjustable
- > Includes fixed mount clamp
- > Offered with Motion Concepts Standard, Stealth All Positioning and Comfort Plus as well as Onyx Auto Style headrests

MOTION CONCEPTS STANDARD HEARDEST



- > Easy to clean self skinning polyurethane foam
- > Gentle lateral contouring
- > 6½" wide x 7" high

STEALTH ALL POSITIONING HEARDEST



- > Soft, comfortable layered foam on coated metal shell
- > Removable, washable cover
- > Four pad sizes: Large and X-Large

NINO - STEALTH COMFORT PLUS HEARDEST



- > Soft, comfortable layered foam on coated metal shell
- > Removable, washable cover
- > 6" wide x 3¾" high

STEALTH COMFORT PLUS HEARDEST



- > Soft, comfortable layered foam on coated metal shell
- > Removable, washable cover
- > Four pad sizes: 8", 10" or 14" wide x 6" high

ONYX AUTO STYLE HEARDEST



- > Onyx leatherette cover
- > Auto-style contouring
- > 8" wide x 8" high x 4" thick