

You say you want an evolution ...

Well, you know, we all want to change the world. And we're doing what we can — one client at a time. Thanks to the new UltraLow CG Tilt system — which evolved directly from our highly regarded Low-Pro and is designed especially for the latest compact power bases — you get the unprecedented flexibility and adjustability you need to provide a customized, perfect fit for each client!

The innovative UltraLow features:


- ▶ **Whole system adjusts on base** — new, slotted channels provide easy fore and aft adjustability (3" - 5").
- ▶ **Independent, telescoping front rigging** and sliding seat pan for easy seat depth adjustment.
- ▶ Sleek design allows use of **maximum battery size** that power base can accommodate.
- ▶ **Extended CG shift** (up to 7") for weight distribution in bariatric or vent applications.
- ▶ **Seat-to-floor heights** as low as 17.25" (with tilt and recline).
- ▶ **Fully modular** — CG Tilt (55°) and/or Recline (174°) and/or Power Elevating Seat systems (up to 7").
- ▶ **Easily interchangeable** from base to base (see list below*).

So get in on the evolution. Check out Motion Concepts' new Ultra Low power positioning system today and give your clients more adjustability, versatility and lower STF's than they ever imagined on a CG Tilt. While it may not actually change the world, it can make your client's world a lot more comfortable.

UltraLow modular CG Tilt



Pictured: UltraLow on Freestyle F11 Base

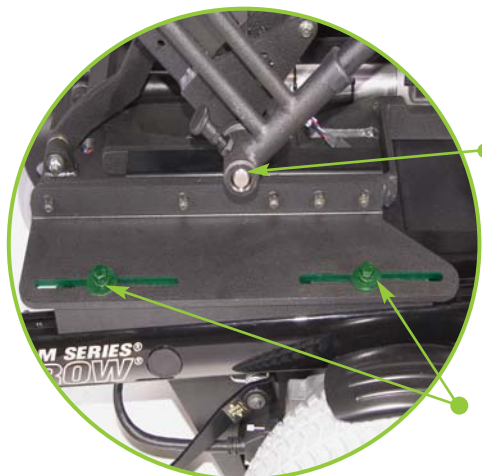

Now up to 7" of CG shift.

Don't lift a thing — the entire system provides simple and fore and aft adjustability with just 1 tool.

*Available on these bases:

F11	Q600
M11	Q610
P220/P222	Q6000
	1103
	1107
Storm	1113 (ATS)
TDX	Dynamo (ATS)
M71	1122
M91	1402

Concepts Motion

US: 700 Ensminger Rd., # 112
Tonawanda, NY 14150
Toll Free Tel: 1.888.433.6818
Toll Free Fax: 1.888.433.6834

Canada: 84 Citation Dr.,
Concord, ON L4K 3C1
Toll Free Tel: 1.866.748.7943
Tel: 905.695.0134
Fax: 905.695.0138

www.motionconcepts.com