

www.motionconcepts.com
U.S.
700 Ensminger Rd.,
Suite 112
Tonawanda, NY 14150
TEL 1.888.433.6818
FAX 1.888.433.6834


Canada
84 Citation Dr.,
Concord, ON
L4K 3C1
TEL 905.695.0134
FAX 905.695.0138

Highly sensitive molded visco-elastic foam top layer.

Foam sweet foam ...

Everyone knows there's no place like home for that good old fashioned feeling of comfort and security. But for many people with disabilities, their wheelchairs sometimes feel like their 2nd homes. And those are just the people Motion Concepts created the Matrix-V cushion for.

The new Matrix-V cushion is made of foam. That's right, foam. Because for true, down-home, comfort that's sweet as can be, there's really nothing like foam. The Matrix-V is an anatomically shaped, foam cushion that provides effective and reliable pressure reduction along with the skin protection and positioning clients need — without any of the hassles of maintenance or adjustment. It's no wonder the Matrix-V has quickly become the cushion of choice for many dealers and therapists.

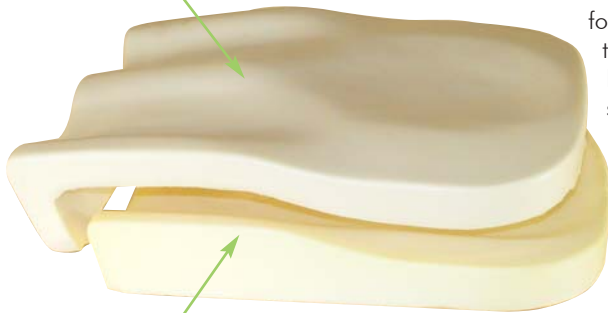

Anatomically contoured, lightweight, molded HR polyurethane foam base.

Best seat in the house

The Matrix-V cushion is created from layered HR and highly sensitive, visco-elastic foam. This amazing combination provides an excellent overall surface contact area (pressure = force/area; pressure mapping study available on request). Plus, clients get great pressure re-distribution as the shaping of the foam increases weight bearing on the hips and thighs, relative to the ischials and sacrum. So, they have the secure feeling of stability, along with skin protection.

There's also a pre-ischial contour which provides excellent positioning — it helps prevent sliding, reduces shear and helps clients maintain an upright position. The Matrix-V also provides outstanding stability through its trochanteric shelf, which helps keep hips level and prevent pelvic obliquity.

The Matrix-V cushion also keeps users cool and dry with a special environmental cover that wicks moisture away. It also comes ready to go! There's no special maintenance needed, and clients get a full 2 year warranty — although it usually lasts *much* longer.

So, check out the lightweight Matrix-V cushion from Motion Concepts today. It's available for quick delivery in standard sizes from 14" to 22" wide, and in 16", 18", 20" standard seat depths. The Medicare billing code is E0192, and facility demos are available by request.

Clearly, for effective and reliable pressure reduction seating, there's really nothing sweeter than foam.

Specially designed three-layer outer cover.



Durable polyurethane coating seals and protects foam