


✓ **MaTRx Flo-tech cushion & MaTRx Flo-tech Fit Kit**

Flo-tech Cushion & Fit Kit

Products shaped around people

MaTRx Flo-tech is a market leading fluid pressure reduction cushion designed especially for users with highest risk of skin breakdown.

- ▶ Sizes from 14" - 20" Wide x 14" - 21" Deep
- ▶ Weight capacity 300 lb
- ▶ HCPCS Code: K0734
- ▶ ADP Code: SES1M0030

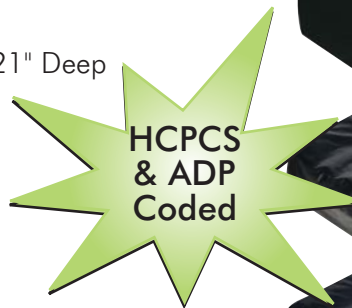

matrx
flo-tech




Stretchy, moisture resistant, breathable outer cover – includes hook and loop fasteners, lifting strap and side pocket

Moisture resistant "Platilon" inner liner protects foam from moisture

Unique dual layer fluid sacs provide immersion to sensitive bony prominences and prevent "bottoming out"

Comfortable, molded foam base includes pre-ischial contour and trochanteric shelf

Flo-tech Fit Kit Components

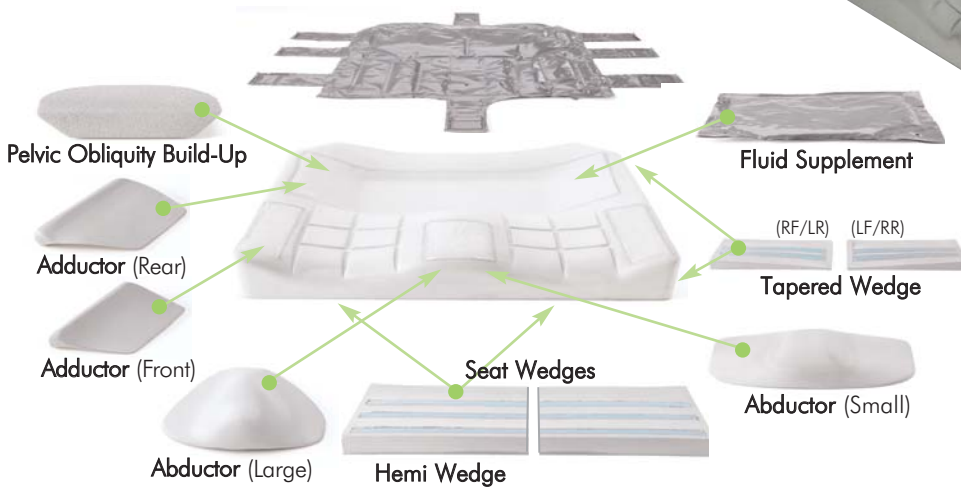

Seat Wedge – Increases degree of pre-ischial contour angle to help prevent sliding

Hemi Wedge – May be used to provide increased support to one lower extremity - useful for clients who use foot propulsion for mobility or with hip flexion discrepancies

Tapered Wedge – Accommodates mild pelvic rotation, mild hip flexion discrepancies or increases adductor contour

Adductor (Front) – Provides increased adductor contour to thighs

Adductor (Rear) – Provides increased adductor contour to hips

Abductor (Large) – Provides 2" increased height of abductor contour

Abductor (Small) – Provides 3/4" increased height of abductor contour

Pelvic Obliquity – Accommodates obliquities up to 1" - for correctable obliquities use under low side of pelvis and for fixed obliquities use under high side of pelvis

Fluid Supplement – May be placed laterally to accommodate mild pelvic obliquities, sacral placement for additional sacral support/pressure reduction, bilaterally for increased ischial immersion.



US: 700 Ensminger Rd. Suite 112
Tonawanda, NY 14150
Toll Free Tel: 1.888.433.6818
Toll Free Fax: 1.888.433.6834

Canada: 84 Citation Dr.,
Concord, ON L4K 3C1
Toll Free Tel: 1.866.748.7943
Tel: 905.695.0134
Fax: 905.695.0138

www.motionconcepts.com