

INVACARE® MATRX®

LIBRA™ CUSHION



The **Invacare® Matrx® Libra™ Cushion** is designed with the highest level of skin protection, positioning and adjustability for unique user needs. It is specially contoured to provide balance and stability to the pelvis, optimizing body function. Recommended for all users, including those at highest risk of skin breakdown.

Call to order

1.800.333.6900
invacare.com

1.888.433.6818
motionconcepts.com



matrx®
SEATING SERIES

Outer Cover – Moisture-resistant, high-stretch, breathable Startex™. Hook and loop fastener strips, lifting strap, convenient zippered side pocket with easy-pull zipper tab. Optional Spacetex 4000U™ cover available for extra breathability, or Infection Control cover made with wipeable Reverse Startex.

Inner Liner – Highly flexible Platilon liner provides an extra layer of moisture protection.

Fluid Sac – Dual-layer, sealed compartments for maximum protection around bony prominences. Non-temperature sensitive, no kneading or maintenance required.

Foam Base – Anatomically contoured to optimize pelvic support and positioning. High resilience foam made with Ultra-Fresh™ anti-microbial and odor protection.


Ab/adductor contour promotes comfortable midline leg position.

Pre-ischial contour prevents sliding to reduce shear.

Trochanteric shelf promotes level pelvis, lateral stability.

Extended rear cushion radius helps prevent posterior pelvic tilt.

Waterfall front edge accommodates tight hamstrings, foot propulsion.

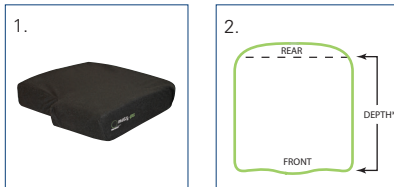


1. Modification

Leg length discrepancy cuts and depth reductions available, outer cover tailored to fit.

2. Actual cushion depth is 1 1/2" longer than stated depth due to rear cushion radius.


HCPCS Code: E2624/E2625

Libra Cushion Fit Kit Accessories: Designed to enhance immersion, correct or accommodate orthopedic issues.

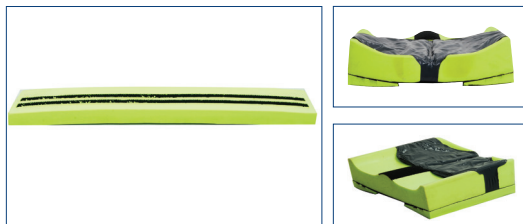


Lateral Wedges

Creates extra depth of contour for enhanced ischial relief

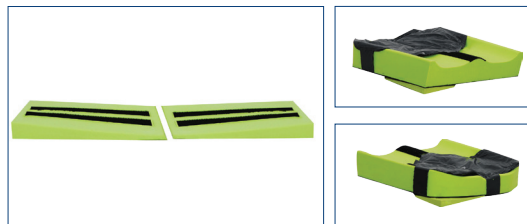


Tapered Wedges

Enhances adductor support, addresses pelvic rotation/obliquity

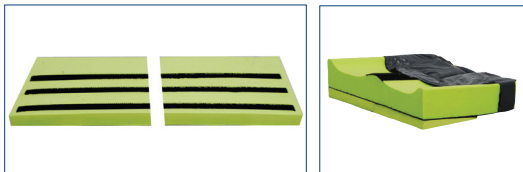

Anterior Seat Wedges

Enhances pre-ischial contour to prevent sliding

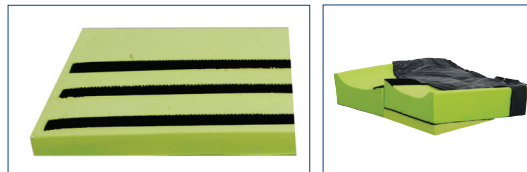


Hemi Wedge

Accommodates hip-flexion limitations, foot propulsion



Fluid Supplement

Unique strap attachment for targeted application

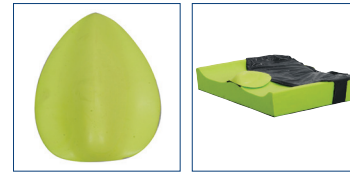

Abductor

Promotes midline leg position

Invacare Corporation
www.invacare.com

USA
One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada
570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

Motion Concepts
www.motionconcepts.com

USA
700 Ensminger Rd., Suite 112
Tonawanda, NY 14150
(888) 433-6818
Fax: (888) 433-6834

Canada
84 Citation Dr., Unit 1-7
Concord, ON L4K 3C1
(866) 748-7943
(905)-695-0134
Fax: 905-695-0138

© 2016
Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols™ and ©. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Matrix is a registered trademark of Motion Concepts. Startex is a registered trademark of Star-Tex Corporation. Spacetex 4000U is a trademark of FR Systems International. Ultra-Fresh is a trademark of Thomson Research Associates. Specifications are subject to change without notification.
Form No. 14-399
Rev. 09/16
160957



Model #	Width	Depth	Height at Abductor/Adductor	Height at Leg Trough	Product Weight	Weight Capacity	Limited Warranty
LC (Standard)	14" - 20"	12" - 22"	3.35"	2.35"	4.0 lb. and up	300 lb.	3 Year
LC (Heavy-Duty)	21" - 24"	18" - 22"	3.60"	2.60"	8.0 lb. and up	500 lb.	3 Year

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program and should consult an attorney or other advisor to discuss specific situations in further detail.

