

# Unbeatable



## INVACARE® AVIVA® STORM RX™ POWER WHEELCHAIR





## INVACARE® AVIVA® STORM RX™

# Our power wheelchairs **take comfort to a whole new level**

Our **AVIVA STORM RX** power wheelchairs are designed to combat the challenges of everyday modern life. At home, at work, out and about, you'll be ready to face what the world throws at you. It's compact, adjustable (so it'll work with you, not against you) and is our most comfortable ride yet.

### **Why choose Rear Wheel Drive?**

- ▶ **Glide over aggressive terrain** with precise control especially at higher speeds due to the optimized RWD suspension system.
- ▶ **Confidently maneuver around your environment** - with only anti-tips behind you, you can handle turns and tight areas with ease.
- ▶ **Unique flared front ends** let you tuck your feet under your chair without any obstructions.



# COMFORT



# FREEEDOM



AVIVA STORM RX base

**GO**

wherever  
you want

You'll power up inclines and cover bumpy terrain like they're smooth pavement and keep up with the pace of modern life. Introducing the new **Invacare® AVIVA® STORM RX™ Power Wheelchair**. Easy to set up, easy to service, and easy to love. It's easy to see why it's Unbeatable.

## Base Technology

### ► Suspension

It gives you great control, precise traction and extra comfort. When the anti-tippers are engaged going over bumps, they push the drive wheel down into the ground to aid traction. It's a compact chair that can really go the distance.

### ► LINX® G-TRAC® technology

The digital gyroscopic technology detects and corrects even the smallest deviation from your intended path, so you'll go exactly where you want to go without having to tweak the control direction. It's especially handy if you use any alternative control inputs with your chair.



# SMOOTH RIDE





AVIVA STORM RX provide  
the perfect balance

# Stable on uneven ground and can tackle almost any terrain

Because all four wheels always stay in contact with the ground, you'll enjoy a smooth and comfortable drive, even when driving on aggressive terrain.

## Base technology

### ▶ Up hill

Going up hills won't be a struggle - you'll power up them, safely and effortlessly due to the DuraWatt™ motors and the way the suspension puts down the power.

### ▶ Silent travels

Quieter ride in the base thanks to the shocks and elastomer dampeners in the caster bearing housing and anti-tippers.

**To truly appreciate the quality of the AVIVA STORM RX, you'll have to take a test drive and experience it for yourself.**

# Slim line design

for those tight spaces

**24"**      **16.75"**

wheel to wheel width      inches Seat to floor height\*

## Compact turning circle

The petite base on the AVIVA STORM RX has a compact turning radius. The flared front end and tuck under capability of the center mount front riggings contribute to the compactness of the base.

Turn in small spaces such as narrow hallways and shopping aisles with ease.

- ▶ Home
- ▶ Local coffee shop
- ▶ Shopping center elevators
- ▶ Office hallways

Navigate through confined areas with ease, while feeling safe and secure in your powerchair.

\*Invacare® AVIVA® STORM RX Power Wheelchair with Ultra Low Maxx







←—————→

## No compromise on distance

Smaller bases used to be hampered by the distance you could travel on them, not anymore – you can have both.





 **STORM** RX20





# It's tailor made for you

Pick from a range of 8 high gloss colors for the contoured shrouds and 11 high gloss colors for the casters and drive wheel rim inserts to match your personality.

Our power wheelchairs have sleek lines and a modern design, so you can feel confident being seen out and about in your chair.



### Power Wheelchair Guide

Use our Power Evaluation App to play with color combinations and visualize your chair before you order.

### Available colors



Obsidian Black Ruby Red



Peacock Blue Purple Passion



Lime Twist Morning Mist



Orange Marmalade Pink Flamingo



Banana Split Arctic Blue



Denim Navy



# LIGHTTTT



# When nights set in

Lights help illuminate the pathway in front of you with our new standard lighting system.

With a bright LED built right into the front caster fork housing and a rear LED strip built into the rear shroud.

# THE NEW WAY


Seat to floor height

# Get closer

**The AVIVA STORM RX has an impressively low seat to floor height so you can get even closer to tables and chairs while using 14" drive wheels and Group 24 batteries.**

It also makes it easier for you to access wheelchair accessible vehicles (WAV).

**16.75"**

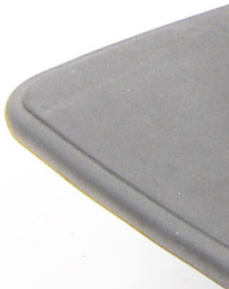







Reach unbeatable heights

# Elevate feature

allows you to  
rise up to 12  
inches

**Whether chatting with friends at a natural eye-level or reaching for objects on shelves.**

▶ **Discreet and dynamic**

Thanks to a sleek elevation mechanism, you can alter the height of your system discreetly, whenever you want thanks to the quiet speed adjustable actuators.

▶ **Moves with you**

There's no need to stop your chair prior to elevating or descending; the simply smart AVIVA STORM RX allows you to utilize the elevate feature while you are driving.





# RISE UP

POWER POSITIONING

# Ultra Low Maxx



We've paired the AVIVA STORM RX's superior base with the Ultra Low Maxx - a high end seating system that makes sure you enjoy perfect posture while feeling comfortable in your chair.

▶ **Natural pivot points**

Natural pivot points is the industry speak for our Ultra Low Maxx feature. It means your head, shoulders, hips and feet stay in contact with the chair, mimicking your body's natural pivot points. It's a snug fit that will leave you feeling safe in the chair.

▶ **Supported at all times**

No matter what dynamic positioning movements the seating system makes, you will remain comfortably supported at all times as the Ultra Low Maxx moves precisely with you.

▶ **Built-in rail**

This allows for accessories to be added to your powerchair depending on your needs and wants.

▶ **Pick your chair**

- Width: 14-22 inches
- Depth: 14-23 inches
- Backrest height: 18-27 inches
- Armrest height: 9-18.75 inches



50° 170°

tilt

recline



# SMART CONTROL

# 3.5"

color touch screen







## Built-in Bluetooth

Connect to your smart phone or laptop.

## Direct Access

Jump straight into the menus - you have direct access to the profiles and functions

## Touchscreen technology

Swipe or tap our color touch screen and be in control of where you go with our REM400 control

## Ambidextrous

The controls are adaptable for both left and right-handed use. If a standard joystick doesn't work for you, you can customize your chair with a specialty control

## Clear speed gauge

On the top left of the control screen is a speedometer, indicating the speed of your chair as you are driving it. An odometer is also present to record the distance you have travelled

## Audible Beeps

Turn on audible cues to get an audio reminder of what function and profile you are navigating to

Color touchscreen joystick available as an option.

# CONTROL

## Specialty Controls



Product	Description	Type
	<p><b>ASL 110/165 Fusion with Bluetooth</b></p>	<p><b>Head controls</b></p>
	<p><b>ASL 102 Uno single switch scanner with Bluetooth</b></p>	<p><b>Single switch scanner</b></p>
	<p><b>ASL 128 Molecule joystick with Bluetooth</b></p>	<p><b>Joystick</b></p>
	<p><b>ASL 104 Atom with Bluetooth</b></p>	<p><b>Head controls</b></p>

ASL specialty controls are available as off-chair purchases only. Please call ASL directly to order specialty control options.







# Travel out and about with confidence!

> **MyLiNX** 

MyLiNX is an app that puts a stop to the run of the mill concerns you might have about battery life, and brushes aside nagging worries about something going wrong with your chair. It shows you vital information about your chair's drive time, battery life and overall usage.





# PIECE OF MIND





# DECISION TIME





Making your mind up

# Take our chair for a spin

Choosing a powerchair is a big decision and not always an easy one. We're at the end of the phone to help, we'll explain the benefits of each model so you can find the one that's right for you. If getting to a dealer to test out a chair is tricky, let them know - they might be able to visit you instead.

**Now you've seen all the AVIVA STORM RX features, it's time to try it for yourself.**

The only way you'll really know if you've found the right chair for you is to try it out.

You can 'test drive' all our power wheelchairs.

**Arrange your demo at your local dealer today.**

# Unbeatable



Design. Technology. Performance.

**INVACARE.COM/UNBEATABLE**

**Invacare Corporation**  
[www.invacare.com](http://www.invacare.com)

**USA**

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

**Canada**

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

**Motion Concepts**

[www.motionconcepts.com](http://www.motionconcepts.com)

**USA**

700 Ensminger Rd. Suite 112  
Tonawanda, NY 14150  
(888) 433-6818  
Fax: (888) 433-6834

**Canada**

84 Citation Dr., Unit 1-7  
Concord, ON L4K 3C1  
(866) 748-7943  
(905) 695-0134  
Fax: (905) 695-0138

©2021

Invacare Corporation.  
All rights reserved.  
Trademarks are identified  
by the symbols ™ and  
®. All trademarks are  
owned by or licensed  
to Invacare Corporation  
unless otherwise noted.  
Specifications are subject to  
change without notification.  
Form No. 20-200 201009

