

SMART Just GOT SMARTER!

Introducing the Motion Concepts' **SMART SHIP** program for Invacare TDX SP2!

Effective June 2018, order your Motion Concepts UltraLow Maxx power positioning system with an TDX SP2 base — directly through us in one simple “smart” step.

Here's what our **SMART SHIP** program means to you:

- ✓ *Invacare's TDX SP2 base is available directly from Motion Concepts*
- ✓ **ONE CALL DOES IT ALL . . .**
ONE order form, ONE purchase order and ONE invoice!
- ✓ *Let our seasoned customer service rehab specialists help you navigate the order form for BOTH products*
- ✓ *Installed in our facility by our highly trained technical specialists*
- ✓ *Motion Concepts' knowledgeable technical support for BOTH products!*

Access . . .

our full line of innovative Power Positioning Products

and ALL the modularity, functionality and durability you have come to expect

from Motion Concepts!!

Call to order

888.433.6818
motionconcepts.com



OUTBACK CANTILEVER RECLINE ARM



SMART Just GOT SMARTER!




© 2018 Motion Concepts, LP. All rights reserved. Trademarks are identified by the symbols ™ and ©. All trademarks are owned by or licensed to Invacare Corporation or its subsidiaries unless otherwise noted. Matrix is a registered trademark of Motion Concepts. Specifications are subject to change without notification. Form No. 17-355D Rev. 06/18

SINGLE POST BACK



ENDOMORPH BACK CANES

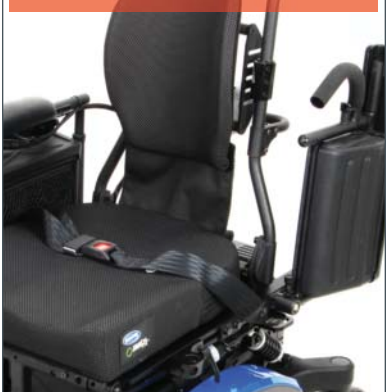

O2 HOLDERS & VENT TRAYS




M290-ATTENDANT CONTROL WITH COLOR TOUCH SCREEN




POWER DOWN TRANSFER ARM



FULL LINE OF INVACARE® MATRX® SEATING

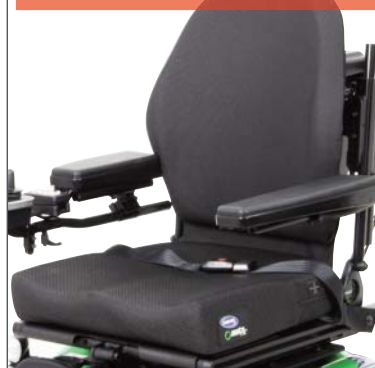

45° POWER ANTERIOR TILT

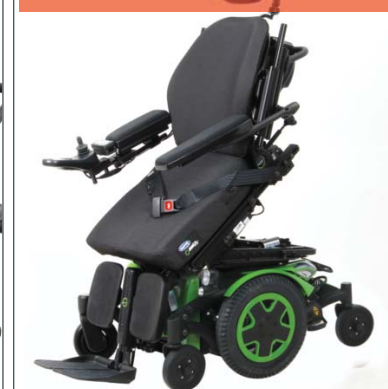

OUTBACK, CANTILEVER RECLINE ARM
The NEW, highly-adjustable and ergonomically designed low-profile OutBack arm offers streamlined design and ability to accommodate most complex arm positions. Arm has width (up to 4" range per side: 3" in and 1" out) adjustment as well as depth and angle adjustment. The OutBack feature means that the arm will always flip back for side transfers, no matter what the position is.

ENDOMORPH BACK CANES
Endomorph back canes are designed to suit clients who are wider at their hips than at their trunk. Available in 1", 2" or 4" seat variances (6" total with arm spacers) and offers clearance for redundant tissue at hips/lower back.

M290 - COLOR, TOUCH SCREEN ATTENDANT CONTROL & PROGRAMMER
Easy to follow color, touch screen attendant control and programmer that controls all power positioning functions:

- Hand-held attendant control of up to six power positioning functions
- Easily set individual actuator speeds in both forward and reverse directions
- Configure installed mechanical switch in either toggle or direct mode, no more external harness or complicated programming
- Program limits and lockout angles through M290
- No need for individual hand-held programmer, everything you need to do can be accomplished through M290 for seating.

POWER DOWN TRANSFER ARM
The weight bearing (up to 200 lbs) transfer arm is available in a power or manual version, allowing clients to slide/shuffle along seat pan onto armrest for convenient transfers.

INVACARE® MATRX® SEATING
The Invacare Matrix Seating line provides award-winning seating solutions from cushions and backs to headrests and other accessories for a wide range of client needs.

O2 HOLDERS & VENT TRAYS
Motion Concepts® offers Articulating & Non-articulating Vent Trays, IV Poles and O2 holders in various configurations, designed to meet the individual client needs.

SINGLE POST BACK
An attendant operated, manual release lever provides back angle adjustments from 73° to 109° in 6° lockable increments. Back folds forward for easy transportation. Offered with our standard steel Rehab Back, optional Contour or High Backs or with an interface plate for use with any of our Matrix backrests.

45° ANTERIOR TILT
A separate power module designed to allow up to 45° of Anterior Tilt, while maintaining 45° of posterior tilt capabilities. This module is intended for use with Ultra Low Maxx Power Positioning system, and is effective in day-to-day activities.

Call to order
888.433.6818
motionconcepts.com

