

MATRIX®
E2™ BACK
SERIES



The NEW Invacare® Matrix® E2™ Back Series has taken back support to a new level! Introducing the new E2 Lateral back featuring our popular BackBone PSIS/Sacro lumbar support. Leading edge Ripple Foam technology is designed to provide temperature control and enhanced immersion for our most comfortable support ever.


Call to order

US: 1.888.433.6818 | CAN: 1.866.748.7943
motionconcepts.com



matrix®
SEATING SERIES

RIPPLE FOAM

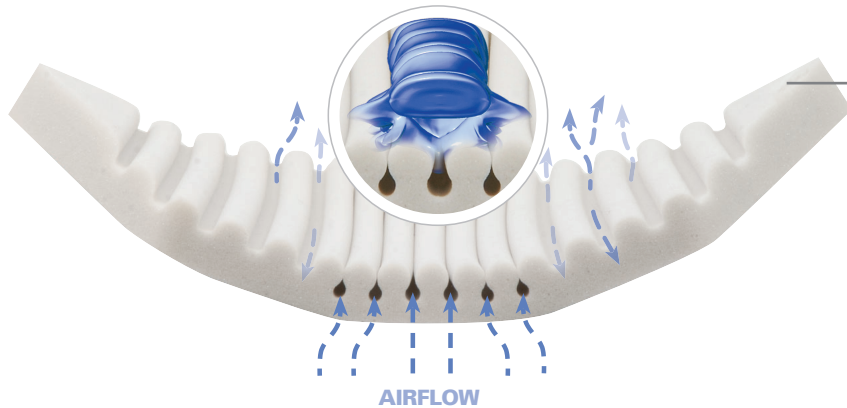

- Airflow channels provide temperature and moisture control
- Enhanced immersion and pressure redistribution
- Extra protection and shear reduction in spinal relief area
- Allows perfect fit of foam to back contour

HARDWARE OPTIONS

Mounting Pin Interface: Height and Width adjustable

Adjustments: Angle (40°), depth (2" in front or behind back canes), can be offset for rotation

Fits: 1" back canes, inserts for 3/4" & 7/8" back canes included

MiniSet
Standard on
6" - 16" Length
E2 Series backs



EasySet
Standard on
18" & 20" Length
E2 Series backs



OPTIONAL

6" - 16" Length E2 Series backs



1. Fixed MiniSet


- Low profile design compatible with rigid chairs, recommended for tilt/recline applications
- Requires tools for removal

2. Heavy-Duty MiniSet

- Recommended for high tone or heavy use applications
- Quick release

OPTIONAL

18" & 20" Length E2 Series backs

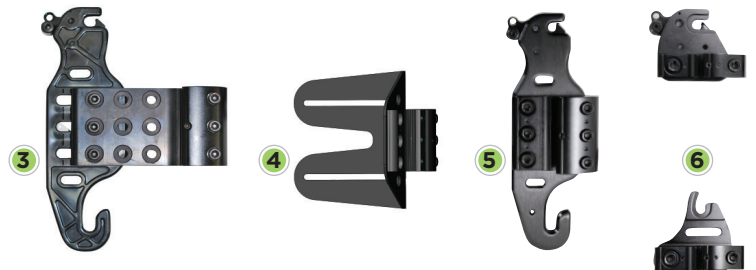

3. Extended Range EasySet

- Compatible with E2S/E2LB/E2D backs lengths 12" to 20" long
- Provides additional 2" depth adjustment

4. Fixed EasySet

- Recommended for tilt/recline applications
- Requires tools for removal

5. Heavy-Duty EasySet

- Recommended for users with high tone or heavy intended use
- Quick release

6. Split

- Allows clearance of frame obstacles, quick release
- Compatible with EasySet Pin Interface only

E2S

Standard Contour Back



E2615




new E2LB

Lateral Back



E2620

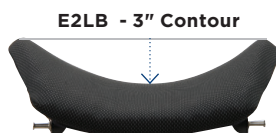

E2D

Deep Contour Back




E2620




Easy slide mounting pin
Assembly




Back model width 1" >
wheelchair back frame width

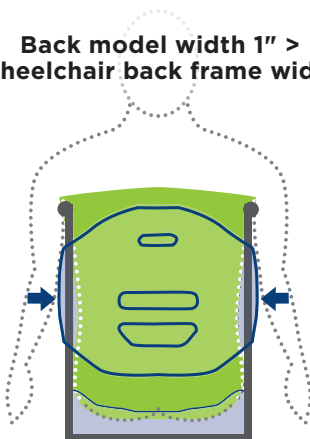

Back model width =
wheelchair back frame width

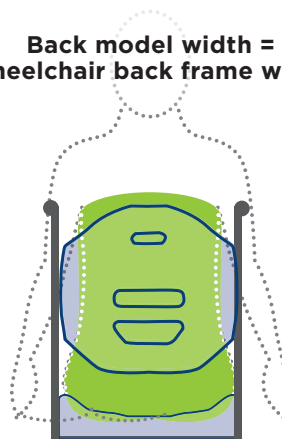

Back model width 1" < to 2" >
wheelchair back frame width

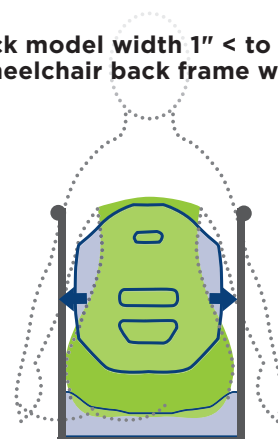

LENGTH

WIDTH

	14"	15"	16"	17"	18"	19"	20"
FITS W/C	13"-16"	14"-17"	15"-18"	16"-19"	17"-20"	18"-21"	19"-22"
6"	●	●	●	●	●	●	●
8"	●	●	●	●	●	●	●
10"	●	●	●	●	●	●	●
12"	●	●	●	●	●	●	●
14"	●	●	●	●	●	●	●
16"	●	●	●	●	●	●	●
18"	●	●	●	●	●	●	●
20"	●	●	●	●	●	●	●

• All E2 back models available in sizes listed

• Up to 300 lb weight capacity

● 6", 8" and 10" length backs include MiniSet hardware

● 12", 14" and 16" length backs include MiniSet hardware – EasySet hardware may be substituted at no charge

● 18" and 20" length backs include EasySet hardware

See left hand page for additional hardware options

ACCESSORIES




Inside Mount

Outside Mount



EasyFix Laterals

Height adjustable with user seated in back.

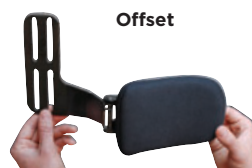

Offset

Standard

Swing Away Laterals

Depth Adjustable pad. Available in 4 pad sizes.



BackBone

Contoured PSIS/Sacrolumbar Support.



PostureFlex

Comfortable, quick release abdominal support. Available in 2 sizes.



Quad Release Tab

Quad release tab for PostureFlex.



Protect Slipcovers

Infection control, washable Startex slip cover.



Matrix Privacy Flaps

Easy attachment to all Matrix Backs with existing sewn on hook & loop. Available in 4 sizes.



Chest Harness Interface

Allows precise placement of Shoulder Strap Guides or Chest harnesses. Available in 2 sizes.



E2HDI

Headrest Interface

Standard interface plate included with all E2 models 14" and up lengths. E2HDI heavy duty interface plate available separately - recommended for extreme tone or heavy use applications, or mounting non-Matrix headrests.



Motion Concepts
www.motionconcepts.com
info@motionconcepts.com

USA
700 Ensminger Rd., Suite 112
Tonawanda, NY 14150
(888) 433-6818
Fax: (888) 433-6834

Canada
84 Citation Dr., Unit 1
Concord, ON L4K 3C1
(866) 748-7943
(905)-695-0134
Fax: 905-695-0138

© 2023Matrix is a registered trademark of Motion Concepts. Specifications are subject to change without notification.

TRD0720 Rev A