

INVACARE®  
**PINDOT & MATRX®**  
INTEGRATED SEATING

---



**Best of Both Worlds!**

Custom contoured PinDot® Silhouette® and ContourU® integrated with Invacare® Matrx® back shells for perfect fit and precise orientation within the chair. A great choice for those with more challenging positioning needs.

Call to order  
**1.800.451.3553**  
invacare.com



# IT'S AS EASY AS 1-2-3



## 1 CAPTURE THE SHAPE

*NEW PinDot® CLIP wireless technology sets up in seconds, giving you the freedom to capture the molded shape untethered to a laptop. CLIP views in color and solid mesh and scans the mold with digital accuracy for the clearest of detail. With the push of a button, the scan information and your product specifications are emailed to PinDot.*

## 2 ORDER FROM PINDOT



## 3 ONE DELIVERY

*PinDot + Matrix complete*

*PinDot + Matrix + Freedom Designs chair complete and installed\**

*PinDot + Matrix + Rovi® X3 Power Wheelchair and Invacare® TDX® SP2 Power Wheelchair complete and installed\**

*\*On-chair US only*

# SILHOUETTE



Capture a client's unique shape using state of the art technology to create a personalized lightweight foam back to accommodate postural asymmetries and moderate orthopedic deformities. Silhouette cushions are designed to provide comfort, stability and skin protection.

**Foam:** Precise, custom-carved polyurethane foam available in soft, firm or extra firm. A moisture-proof coating protects the foam for easy wipe off cleaning.

**Back Shell Options:** Invacare® Matrix® Elite™, Elite® TR and Elite® Heavy-Duty, Posture Back and Posture Back Heavy-Duty.

**Cover:** Custom designed to match each unique cushion shape and fit securely on the Matrix shell.


**Hardware:** Matrix Hardware allows for simple height, depth and angle adjustments while offering easy removal for portability. Available in fixed or quick release.



Matrix Hardware



Matrix Back Shell

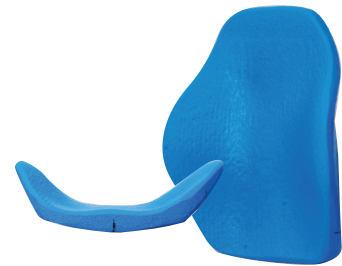

Silhouette Custom Foam Insert/Cover

# CONTOURU & CONTOURU 4U



*Internal ABS Plastic or Aluminum custom contoured lateral reinforcements provide full trunk support like no other cushion*

PinDot ContourU is designed to provide the highest level of support and stability for users with the most complex positioning needs. Multiple options are available to fine tune for extra skin protection and positioning needs.

A moisture-proof coating seals and protects the foam from moisture and allows the cushions to be easily wiped clean. For additional foam protection, add an optional moisture-proof and anti-microbial inner cover.

ContourU 4U is a pre-configured solution for ease of ordering. Simply send your 3D CLIP image, pick your mounting, width, everything else is pre-elected for you. Quicker quotes, orders and delivery.

**Foam:** Designed for those with complex positioning needs, choice of standard or soft foam custom poured polyurethane provides full contact support, even weight distribution and management of orthopedic deformities and muscle tone.

**Back Shell Options:** Matrix Elite, Elite TR and Elite Heavy-Duty, Posture Back and Posture Back Heavy-Duty.

**Cover:** Custom designed to match each unique cushion shape and fit securely on the Matrix shell.


**Hardware:** Matrix hardware provides ease of height, depth, angle and even rotation adjustment, all with user seated in the back support. Several options available in quick release for portability, or fixed, non-removable.



Matrix Hardware



Matrix Back Shell

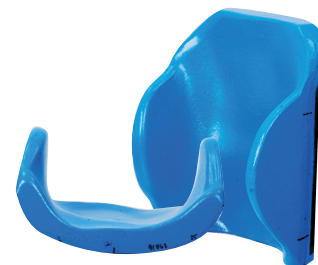

ContourU Custom Foam Insert/Cover

# BACK SHELL OPTIONS



### Elite/Elite Heavy-Duty Back Shell

Lightweight aluminum back shell with scapular cutout and 3" contour depth. Complete with standard and long mounting pins. Standard 12"-20" wide x 10"-20" high (300 lb weight capacity). Heavy-Duty 20"-30" wide x 16"/20" high (600 lb weight capacity).



### Elite TR/Elite TR Heavy-Duty Back Shell

Aluminum back shell including shoulder support, designed for tilt recline applications. Includes standard and long mounting pins. Standard widths 15"-20", heights 14"-20" (300 lb weight capacity). Heavy-Duty 21"-24" wide x 16"-20" high (600 lb weight capacity).



### Posture Back/Posture Back Heavy-Duty Back Shell

Aluminum back shell including shoulder support, designed for tilt recline applications. Includes standard and long mounting pins. Standard widths 15"-20", heights 14"-20" (300 lb weight capacity). Heavy-Duty 21"-24" wide x 16"-20" high (600 lb weight capacity).

# MOUNTING HARDWARE OPTIONS

Easy Set



PBMH

Split



SMH

Heavy-Duty Easy Set



HDMH

Extended Range



ERMH

Fixed Space Saver



FSSH

### Easy Set (PBMH)

Compatible with all Matrix standard backs. 300 lb weight capacity. Adjustable depth 2" in front or behind back canes. Adjustable height, 40° angle, rotation. Quick release portability.

### Split (SMH)

Allows clearance of back cane rigidizer bars or other frame obstacles. Compatible with all Matrix standard and Heavy Duty backs. 500 lb weight capacity. Adjustment range same as Easy Set hardware above.

### Heavy-Duty Easy Set (HDMH)

Required for all Matrix Heavy-Duty backs, recommended on standard backs for users with high tone or heavy use. 500 lb weight capacity. Adjustable depth 2" in front or behind back canes. Adjustable height, 40° angle, rotation. Quick release portability.

### Fixed Space Saver (FSSH)

Recommended for power tilt/recline or tighter width installs. Permanent mount requires tools for removal. Compatible with all Matrix standard and Heavy-Duty backs. 500 lb weight capacity. Adjustment range same as Heavy-Duty hardware above.

### Mounting Pins

**Standard/Long Mounting Pins:** Included with PB and Elite Backs; allows PB Back to fit 1" narrower to 1" wider and Elite backs the same width or 1" wider wheelchairs.

**Optional Extended Mounting Pins (EMP):** Allow interface to wheelchairs up to 2" wider. 300 lb. weight capacity maximum.

Invacare Corporation  
www.invacare.com

### USA

1200 Taylor Street  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

### Motion Concepts

www.motionconcepts.com  
info@motionconcepts.com

### USA

700 Ensminger Rd., Suite 112  
Tonawanda, NY 14150  
(888) 433-6818  
Fax: (888) 433-6834

### Canada

84 Citation Dr., Unit 1-7  
Concord, ON L4K 3C1  
(866) 748-7943  
(905)-695-0134  
Fax: 905-695-0138

### © 2018

Invacare Corporation.  
All rights reserved.  
Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Matrix is a registered trademark of Motion Concepts. Specifications are subject to change without notification.  
Form No. 18-246  
Rev. 08/18  
180878

