


MATRIX®  
**MAC BACK**  
ADJUSTABLE CONTOUR

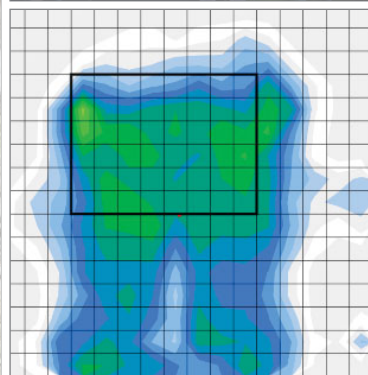

Introducing a whole new approach to bespoke back support, the **Matrix MAC Adjustable Contour Back Series!** Completely configurable for the widest range of body shapes and conditions, and adjustable for changes, the MAC back offers a whole new way to provide custom fit, positioning and comfort!

Call to order

US: 1.888.433.6818 | CAN: 1.866.748.7943  
motionconcepts.com





## LATERAL CONTOUR OPTIONS



**Width Adjustment**

- Allows 1/2" of width adjustment on each side, total 1" growth adjustment



### Height Adjustment

- Lateral contours height adjustable with user seated in the back



### Installation Tool

- Provides ease of setup aligning hardware and lateral contours during installation

### MAC - DUAL LAYER FOAM MODEL



Supersoft/firm HR foam layers for comfort and support



- Supersoft over firm HR double layer foam
- Support and comfort even in tilt/recline applications
- **Recommended for symmetrical postural presentation**

### MAC - HUG MODEL



Soft, conformable StaminaFibre® for maximum immersion and envelopment



- StaminaFibre® insert, machine wash and dry friendly
- Highly immersive micro climate control
- **Recommended for accommodation of mild postural asymmetries**

### SHAPE insert

- Customizable shape interface for accommodation of orthopedic asymmetries



### MAC - HUG MODEL/SHAPE



SHAPE Insert layered under HUG StaminaFibre®



- Positioning interface with configurable HR foam pieces
- Easily adjustable for future changes
- **Recommended for accommodation of moderate postural asymmetries\***

\*(For HUG applications only)

**MAC - CUSTOM: Looking for more customization?** explore the Matrix® MAC Custom Back Support



A bespoke solution, fully configurable to deliver personalized fit, optimal positioning, and lasting comfort.




Custom Order Form


Easy slide mounting pin Assembly




Back model width 1" > wheelchair back frame width

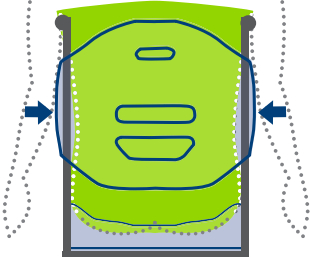

Back model width = wheelchair back frame width

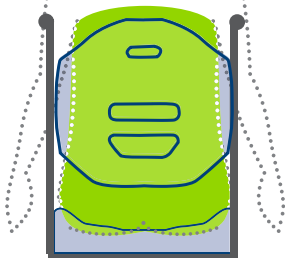

Back model width up to 2" < wheelchair back frame width

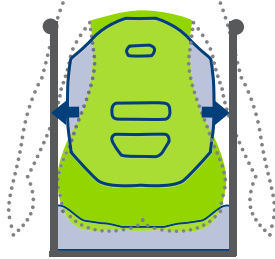

BACK SIZES


WIDTH

LENGTH

|          | 14"     | 15"     | 16"     | 17"     | 18"     | 19"     | 20"     |
|----------|---------|---------|---------|---------|---------|---------|---------|
| FITS W/C | 13"-16" | 14"-17" | 15"-18" | 16"-19" | 17"-20" | 18"-21" | 19"-22" |
| 14"      | •       | •       |         |         |         |         |         |
| 16"      | •       | •       | •       | •       | •       | •       | •       |
| 18"      | •       | •       | •       | •       | •       | •       | •       |
| 20"      | •       | •       | •       | •       | •       | •       | •       |
| 22"      |         |         | •       | •       | •       | •       | •       |
| *24"     |         |         | •       | •       | •       | •       | •       |

\*24" Length is Custom Up to 300 lb weight capacity

OPTIONAL HARDWARE

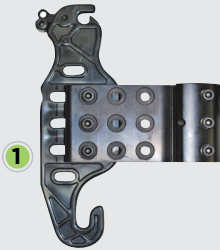


- 1. Extended Range EasySet**
- Provides 4" forward or rear depth adjustment



- 2. Fixed EasySet**
- Recommended for tilt/recline applications
  - Requires tools for removal



- 3. Heavy-Duty EasySet**
- Recommended for users with high tone or heavy intended use

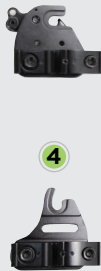

- 4. Split**
- Allows clearance of frame obstacles, quick release



**Fixed Laterals**

Height and width adjustable  
Available in 6 sizes.



**Swing-Away Laterals**

Height and width adjustable  
Available in 4 pad sizes.



**PostureFlex with Tether Strap**

Comfortable, quick release  
anterior trunk (abdominal) supports.  
Available in 2 sizes.



**Quad Release Tab**

Quad release tab for PostureFlex.



**Protect Slip Cover**

Machine-washable, Infection  
Control Slip Cover



**BackBone**

Contoured PSIS/Sacrolumbar  
Support.



**Matrix Privacy Flaps**

Easy attachment to all Matrix  
Backs with existing sewn on hook  
& loop. Available in 4 sizes.



**Chest Harness Interface**

Allows precise placement of  
Shoulder Strap Guides for Chest  
harnesses. Available in 2 sizes.



**Motion Concepts**  
www.motionconcepts.com  
info@motionconcepts.com

**USA**  
700 Ensminger Rd., Suite 112  
Tonawanda, NY 14150  
(888) 433-6818  
Fax: (888) 433-6834

**Canada**  
505 Cityview Blvd. Unit 1  
Vaughan, ON L4H 0L8  
(866) 748-7943  
(905)-695-0134  
Fax: 905-695-0138

© 2025 Matrix is a registered trademark of Motion Concepts. Specifications are subject to change without notification.